



INDIVIDUAL TABLE SOLUTIONS

FOR A MORE ACTIVE
QUALITY OF LIFE!



TRANSFORMING CARE. IMPROVING LIFE



CONTENTS

CHOOSING THE TABLE THAT IS RIGHT FOR YOU

What a difference a table can make

HOW TO GET THE BEST SITTING POSITION

For all people, regardless of mobility

ROPOX VISION TABLES

Functionality and design melted together

ROPOX ERGO TABLES

Guaranteed to match your essential needs

ROPOX GROUP TABLES

For the moments we share

CHOOSING THE TABLE THAT IS RIGHT FOR YOU

When buying a new table, there are certain important features you should consider making sure that the table matches your requirements.

How many will be using the table?

If there is only one user of the table, you should pick a table that best matches the individual user's needs. If there is more than one user of the table, you should choose a table that is easy to adjust in height and a size that fits several users. Furthermore, you may want to consider a table with wheels which you can move around because it will give you more flexibility.

Are there specific needs?

A safe sitting position requires a stable table and chair, especially for people with balancing challenges. If the table has wheels, it is imperative that the wheels have brakes to prevent accidents.

If you have visual challenges or maybe a limited reach, tables with tiltable tops are designed specifically to suit these needs.

Are there sufficient legroom under the table?

Being able to get close to the table top is essential, and there must be good access for a wheelchair user from both sides of the table. All Ropox tables have optimum access possibilities as there are no cross bars underneath the worktop.

WHAT A DIFFERENCE A TABLE CAN MAKE

A good ergonomic sitting position can make a significant difference in your daily life. Children, adults, walking-impaired, people with paralysis and blind people all have different demands, but the flexibility of our tables will provide you with the best possible sitting posture because you can adjust the height, choose from a range of table top sizes and get worktops that tilt.

With the right table, you can be active reading, working or playing and you can get the most out of your daily activities.

TABLE TOP TILT

Tilttable table tops give the possibility to vary your sitting position.



Laptop work – Relaxed back-leaned posture using the front of the table.



Writing – Posture leaning forward using the middle of the table.



Reading – Relaxed posture using the top of the table and bringing it close to you.

TABLE HEIGHT ADJUSTMENT

A height adjustable table makes it easy to adjust the table exactly to your height.



Adult



Person standing

FOR ALL PEOPLE, REGARDLESS OF MOBILITY.

A table is a gathering point in life where we talk, eat, work, play and laugh. It is for all age groups, for walking-impaired, for wheelchair users, for blind people and any other person - disabled or not.

Our tables are created for the moments in life that we share with others, regardless of mobility. That is why you will find our tables in private homes, at schools, nursing homes, rehabilitation centres and hospitals, and it is the reason our products are preferred by the users, their families, carers and therapists.



HOW TO GET THE BEST SITTING POSITION

Staying in the same seated position for more than 1-2 hours will become uncomfortable for most people. The reason is that sitting is one of the worst stresses for your lower back, shoulders and neck and it may cause muscle pain.

The key to a perfect and healthy sitting position is variation. However, if you have a disability, especially if you are using a wheelchair, it can be a challenge to change your sitting posture, but with the right table you will get the needed support to do this.

We have considered all aspects to design tables that allow you to sit comfortably.

- You have extra support for your arms which relieves the stress on your back and shoulders.
- The height adjustment prevents shoulder and neck pain.
- The tilttable worktop brings documents closer towards you.
- Your documents stay in place with the MagRule magnet even when the table top is tilted.
- You are guaranteed to have easy access from both sides of the table as there are no cross bars that limit the available legroom under the tabletop


All the table's details accommodate you in getting the best possible sitting position no matter if you are writing, reading, playing or just unwinding at the table.

VISION TABLE ADVANTAGES

The Ropox Vision series has a very wide table range with a lot of individual adjustments.



Ropox arm supports

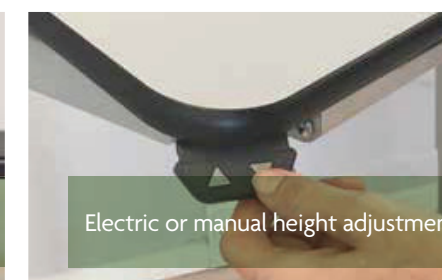

Electric or manual height adjustment



Vision table tops

The Vision table tops have built-in steel plates making you able to place the MagRule wherever you choose.



ROPOX VISION TABLES

FUNCTIONALITY AND DESIGN MELTED TOGETHER

In our opinion, a product must be both functional and beautiful, and the Vision table is the embodiment of such a product. The design is simple and elegant, and it has all the features you need for both work, school, play and training, e.g. adjustable height, a table top that can tilt and a magnet to keep documents fastened to the work top.

With the ergonomic arm supports you can relieve your forearms and they also provide a slight body support. Furthermore, there is plenty of room under the table for a wheelchair user to get close to the table. Many choose the Vision table as their favourite because of the stylish look, user-friendly features and high usability.



The Vision series is ideal for weak sighted due to the almost vertical tiltable table top shortening visual distances.



ROPOX Vision Table

- Table top available with a fixed and tiltable side with the tilt (up to 74°) placed in the middle, left or right side
- Available in two heights (50-70 / 60-90 cm.)
- Three table top sizes (90x60, 120x60, 90x90 cm.)
- Manual or electric height adjustable
- Built-in safety stop on tiltable table top
- No cross bars giving easy access for wheel chair users
- Two built-in guide wheels giving a nice visual appearance



ROPOX Vision High/Low Table

- A table where you can sit or stand
- Electric height adjustable (very large range 62-127cm.)
- Battery solution is available making it cordless
- Tiltable table top
- Steel in the table top - MagRule can be placed arbitrarily to hold documents



The MagRule gives the possibility to place things arbitrarily on the table top due to built in steel plates.



Built-in wheels

Height adjustable

Removable arm supports

Rounded edges

ROPOX VISION TABLE

Max. load capacity: 220 lb
 Height: 23.5" to 35.5"
 Depth: 23.5" or 35.5"
 Width: 35.5" or 47"
 Inner width: 31.5" or 43"
 Frame Colour: Silver Grey


TECHNICAL INFORMATION

VARIANT	A	B	C	D	E
SIZE 35.5 x 23.5					
manual	20-60020-11	2060021-11	2060022-11	2060023-11	2060024-11
electric	20-60020-21	2060021-21	2060022-21	2060023-21	2060024-21

VARIANT	A	B	C	D	E
SIZE 47 x 23.5					
manual	20-60220-11	20-60221-11	20-60222-11	20-60223-11	20-60224-11
electric	20-60220-21	20-60221-21	20-60222-21	20-60223-21	20-60224-21

VARIANT	F	G
SIZE 47 x 23.5		
manual	20-60225-11	20-60226-11
electric	20-60225-21	20-60226-21

COLOUR FOR TABLETOPS:

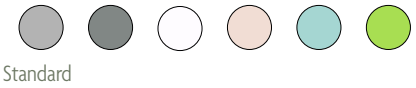

VARIANT	A
SIZE 35.5 x 35.5	
manual	20-60320-11
electric	20-60320-21

ACCESSORIES

Arm support
incl. bag and hook
Item: 20-60520



Arm support lock-able
incl. hook and bag
Item: 20-60517



MagRule
11.8" Item: 20-60545
19.7" Item: 20-60541



Magnet
Item: 50200-000



Castors with brake, without height
increase, table depth 35.5" Item:
20-60589/20-60585



Castors with brake,
increase height by 3"
table depth 23.5" Item:
20-60570/20-60575



Castors with brake, without
height increase,
table depth 23.5"
Item: 20-60573



Castors with brake, without
height increase,
table depth 35.5" Item:
20-60581/20-60579





ROPOX ERGO TABLES



ROPOX Ergo Table

- Manual height adjustable from 56-90 cm.
- Available in two sizes (90x60 and 120x60 cm.)
- Available with full tilt or left/right section tilted with stepped angles up to 54°
- Wide range of accessories



ROPOX Ergo Desk

- Electric height adjustable
- Available with three table top sizes
- Electrical height adjustment from 60-125 cm.

Wide range of accessories



GUARANTEED TO MATCH YOUR ESSENTIAL NEEDS

The Ergo Table consists of four variants that you can adjust manually in height. It has a table top that can tilt up to 54° which is very convenient for weak-sighted and users with limited reach. The supporting ruler will ensure that things will stay in place when the table top is tilted. You can use the arm supports to relieve your forearms and, if required, they can also provide you with a slight body support.

The Ergo Desk is a simple, electric height adjustable table for all.









Whether you need a table for training, teaching or work, the Ergo series would be an option to explore.



ROPOX ERGO TABLE

Max. load capacity:	220 Lbs
Height:	22-35.4"
Depth:	23.5"
Width:	35.5" or 47"
Inner Width:	28" or 41"
Frame Colour:	Silver Grey

TECHNICAL INFORMATION

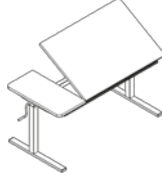VARIANT	A	B	C	D
Size 35.5 x 23.5 Tilt 23.5" / no tilt 11.8"	 20-70010	 20-70011	 20-70012	 20-70013
VARIANT	A	B	C	D
SIZE 47 x 23.5 Tilt 35.5" / no tilt 11.8"	 20-70210	 20-70211	 20-70212	 20-70213

All tables are manually height adjustable from 22 - 35.4"
 An additional height extension can be used for an adjustment interval from 26-39"
 The tabletop can be tilted with stepped angles up to 54° (B to D variant).

ACCESSORIES

Arm support
 incl. bag and hook
 Item: 20-60520



Arm support lock-able
 incl. hook and bag
 Item: 20-60517



Supporting ruler
 Item: 20-70509/35.5"(B)
 Item: 20-70510/47" (B) Item:
 20-70513/35.5"(C+D) Item:
 20-70514/47" (C+D)



Additional height
 extension, 4"
 Item: 20-70503



Document/PC holder
 incl. mounting plate
 Item: 20-70504



Pen holder
 Item: 20-70506



Castors with brake
 without height increase
 Item: 20-70502



4 Brake wheels, increase
 height by 3.5"
 Item: 20-50300



4 Brake wheels, increase
 height by 3"
 Item: 20-50340





ROPOX GROUP TABLES

FOR THE MOMENTS WE SHARE

Meeting with friends, family or colleagues around a table is a common thing in many people's everyday lives, and our range of Group tables are designed to be the natural gathering point whether it is at work, at home, for training or in school.

All Group tables are easy to adjust in height, they are made of quality materials that are easy to clean, and they have plenty of legroom which makes them fit for all users, regardless of the degree of mobility. So, if you are looking for a modern and elegant Group table for your interior which also accommodates for the requirements of disabled users, our range of Group tables will surely be worthwhile looking into.



ROPOX
Vision Group Table

- Manual height adjustable
- Rounded edges and corners
- Smooth surfaces and no dirt traps
- Elegant oval shaped legs



ROPOX
Ergo Group Table

- Manual height adjustable
- No cross bars making it perfect for wheelchair use
- Can be mounted with brake castors




ROPOX
Ergo Multi Table


- Ideal for training situations
- No cross bars
- Different table top cut-out, for individual needs
- Manual height adjustable
- Supplied with brake wheels



Cabinets can be mounted under the 4Single Group table top.

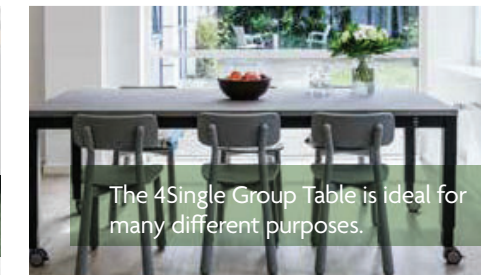

The 4Single Group Table is ideal for many different purposes.

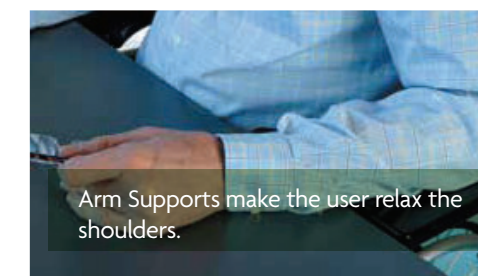

Arm Supports make the user relax the shoulders.

ABOUT ROPOX

FREEDOM AND CONFIDENCE TO BE INDEPENDENT

“People with disabilities deserve an easier life with more self-esteem, quality of life and freedom to safely look after themselves.” This philosophy has been our goal since we set up our business back in 1962. Today, we are more than 60 employees who continue to pursue this mission.

With more than 50 years of experience, we are amongst the world’s leading specialists and producers of assistive devices and furniture for the people with disabilities as well as equipment for ergo- and physiotherapy and for rehabilitation. Being part of the AddLife group, we also have access to a large and strong network of partners who specialise in developing services and products to both the private and public healthcare sector.

Ropo is provider in North American markets by ProCare Innovations .

Find out more on www.procare-innovations.com.



ROPOX TABLES

BEAUTIFUL HEIGHT-ADJUSTABLE TABLES FOR INDIVIDUAL NEEDS



ROPOX KITCHENS

FLEXIBLE AND HEIGHT-ADJUSTABLE KITCHENS



ROPOX BATHROOMS

INNOVATIVE AND EXCITING PRODUCTS IN A BEAUTIFUL COHESION



TRANSFORMING CARE, IMPROVING LIFE