



## THE BATHROOM OF THE FUTURE

- FOR HEALTHCARE AND HOME

FOR A SAFER AND MORE  
INDEPENDENT LIFE



# WASHBASINS

## FOR AN EASIER LIFE

A washbasin must be easy to use, but it should also serve as a support for the user. That is what we keep in mind when we design all our basins: They must be user-friendly and safe.

To make them as user-friendly as possible, they can all be adjusted to meet individual needs, e.g. they are height-adjustable and some can be rotated out in the room so that you can reach them from the toilet if you want to wash your hands while seated.

With our assistive products, you become more independent, and there is of course also the benefit that you will get a stylish bathroom interior that is easy to clean.



Fixed  
Height



Manual



Electric





Solid support grabs for a safe movement

## ROPOX SwingLine



- Primarily suitable for wheelchair users
- The design provides sufficient legroom under the washbasin
- Suitable for sitting users, e.g. when washing hands, brushing teeth etc.
- Its ability to rotate and swing makes the washbasin accessible from all sides
- Can be adjusted in height to meet individual needs
- Cleaning-friendly design without dirt traps
- Dock-in function is recommended to fasten the washbasin

## ROPOX SupportLine



- Especially suitable for walking-impaired
- The integrated, ergonomic grab rails can be used as a support to prevent falling injuries
- Can be height adjusted to meet individual user's needs
- Sideways slidable to minimize the distance when transferring to the washbasin
- The design provides sufficient legroom under the washbasin
- Simple and cleaning-friendly design minimizes bacterial growth

## ROPOX AdaptLine



- Small, elegant, and disability-friendly
- Can be height adjusted to meet individual user's needs
- Highly suitable for wheelchair users with maximum legroom beneath the washbasin
- Easily adapted to existing water and drain connections
- Three models: standard, support, and hospital
- Support Washbasin includes integrated, ergonomic grab rails
- Easy to clean design to allow for high level of hygiene

## ROPOX QuickWash



- Fits into small or adapted bathrooms
- Affordable, manually height-adjustable washbasin
- Simple hi or low position for standing or sitting
- Projects from wall when lower for increased leg room




## TOILET SUPPORTS

### PRODUCTS THAT MAKE YOU MORE INDEPENDENT

Most people want to be independent and handle their daily activities on their own – especially when it comes to personal care and hygiene. However, this may not be possible for people with physical or mental disabilities.

With the right aids and assistive products, you get the support you need, allowing you to utilize the toilet without a helper's assistance. Perhaps, a pair of toilet support arms is all that is needed to get up and reach a standing position. Those with reduced mobility can use the toilet arms to hold on to while sitting down and to assist in getting back into the standing position. The Toilet Support Arms can be adjusted in height and in different angles to suit the individual user's needs and are available in wave or straight versions with varying lengths, and weight capacities. Accessories such as a toilet roll holder and standing/sitting attachment increase the benefits of this assistive product.



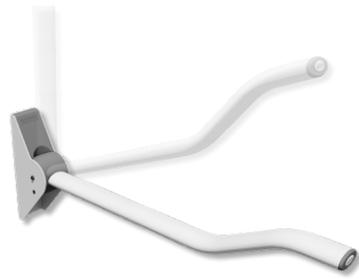

## ROPOX Toilet Support Arms

- Offer support and safety when moving around
- Can be adjusted in various positions
- The oval shape relieves the pressure on your hands and forearms
- Height adjustable and with two integrated heights - optimal support when seated and standing
- Cleaning-friendly – the smooth surface is easy to clean



Provide support at several levels – from a sitting to a standing position.



The standing/sitting attachments provides further support when getting up.



Make it easy to get on and off the toilet.




# TOILET, SINK, SHOWER

## SAFE BATHING AND SHOWERING

For people with limited mobility, the combination of water, soap, and slippery tiles may feel uncomfortable and unsafe, and it may cause an increased risk of falling injuries and accidents. That is why you need solid and stable solutions that make the bathroom safe so that you can feel safe as well.

The solid grabs and shower rails offer a stable support, and if you need to rest, the toilet seat offers great support as well – even in a cramped place. When designing and furnishing your bathroom with the proper aids and assistive products, you can reduce the need for assistance when showering, while gaining the freedom to be more independent.



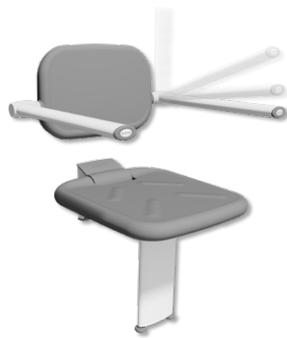

## ROPOX Shower seat

- Reduces the risk of falling accidents
- Height adjustable option
- Non-slip seat with a large surface and a solid grabbing edge
- Hard-wearing, insulating material which does not feel cold on your skin
- Good padding which makes the seat comfortable to sit on and prevents pressure sores
- Armrests and seat can be folded up to a vertical position leaving more space in the room




## ROPOX Shower/changing bed

- Height adjustment to meet the user's individual needs
- The slats have a soft top side, which makes them very comfortable to lie on and a flatter side for easier rolling
- Easy up and down adjustment by using a hand control
- Can be folded up to provide extra space in the room
- Cleaning-friendly as the slats are very easy to remove
- Back support can be mounted



Seated showering  
reduces the risk of  
falling injuries

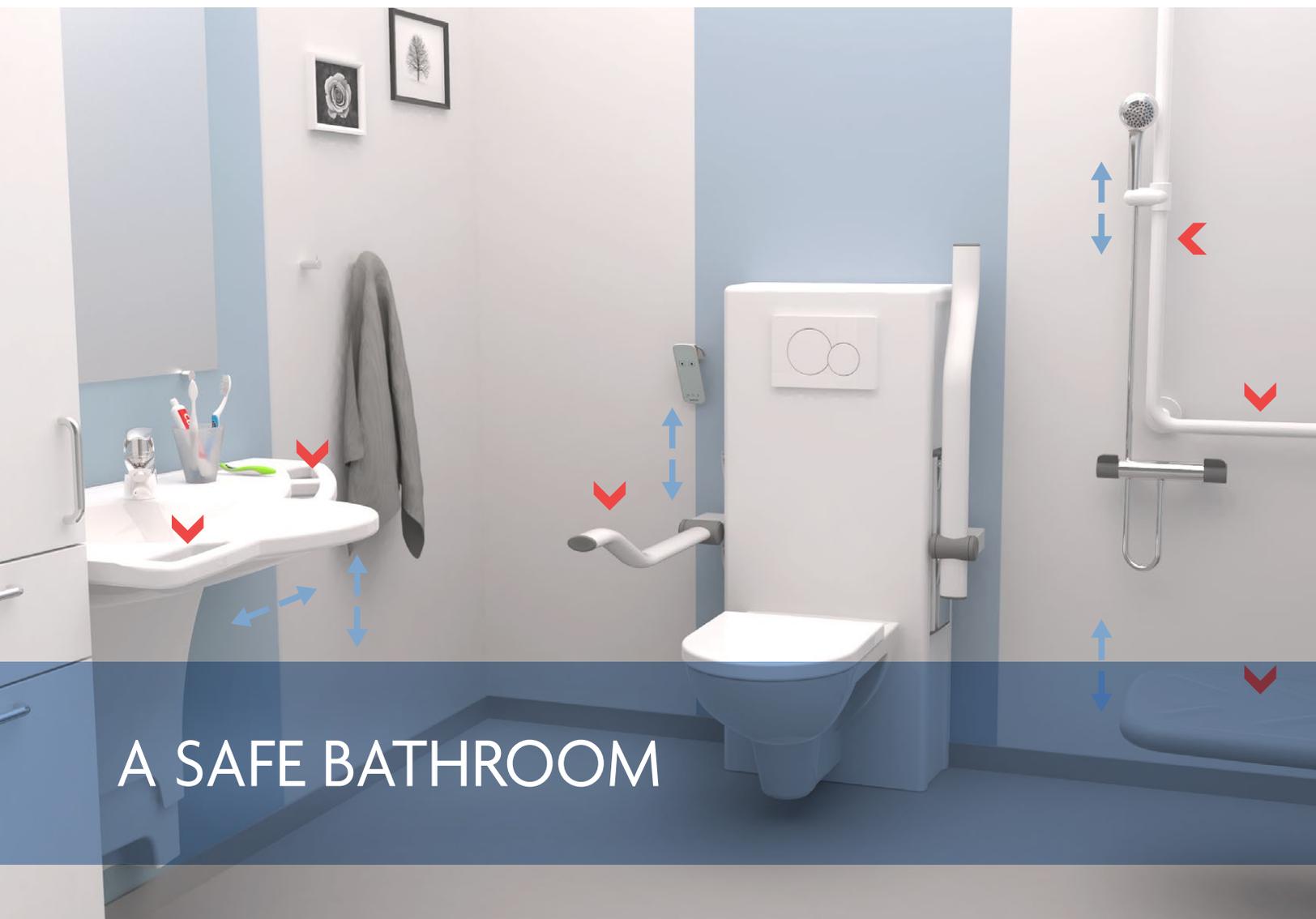

## A SAFE BATHROOM

### FOR THE BATHROOM INTERIOR

Paint an area behind the bathroom furniture and cabinets in a contrasting color (max. 2 different colors), in a luminance contrast of 30% to make it easier for visually impaired to see the white elements.

Choose a flooring material in a light color. A floor with just one solid color makes it easier for people with dementia to orientate themselves in the room because they may perceive dark spots on the floor as holes.

Use anti-slip linoleum on the floor. That makes the room safer - even when the floor is wet.

Choose a sliding door for wheelchair users. It is easier for them to get in and out, it will also provide a larger turning area in the bathroom. Choose a traditional door for people with dementia. They may have difficulties using a sliding door.

Paint the door in another color than the walls. It is easier for visually impaired and people with dementia to identify the door.

■ SUPPORT

■ ADJUST-  
MENT

## FEELING COMFORTABLE AND SUPPORTED

It is very important to furnish your bathroom with many support points and the right aids and assistive products, but other aspects should be considered as well:

- Do you need grab and shower rails in the bathroom for support?
- Do you need a height adjustable toilet and washbasin?
- Are you able to shower in a standing position, or do you need a shower seat?
- Do you need toilet supports at the toilet?
- Do you need more assistance than toilet supports to sit down on and get up from the toilet?
- Would it be an advantage if you could wash your hands and complete personal care while seated on the toilet eliminating the need to move to the washbasin?
- Is the distance between the products minimized to reduce the number of lifts and movements?
- How do you ensure optimal space for a wheelchair user?


When the bathroom is furnished with flexible and user-friendly aids, you will feel safe, and when you feel safe, you are often able to be more independent. For many people, that is the essence of quality of life, and we are very proud to contribute to that.

In our bathroom concept video you will find more inspiration – please find it on [ropox.com](http://ropox.com) under the menu Bathroom concepts.




# BATHROOM DESIGN IDEAS


Designed according to The Danish Building and Property Agency's guidelines for accessible design 2014, Building Regulations BR10 and BR08 and ISO 21542:2012.

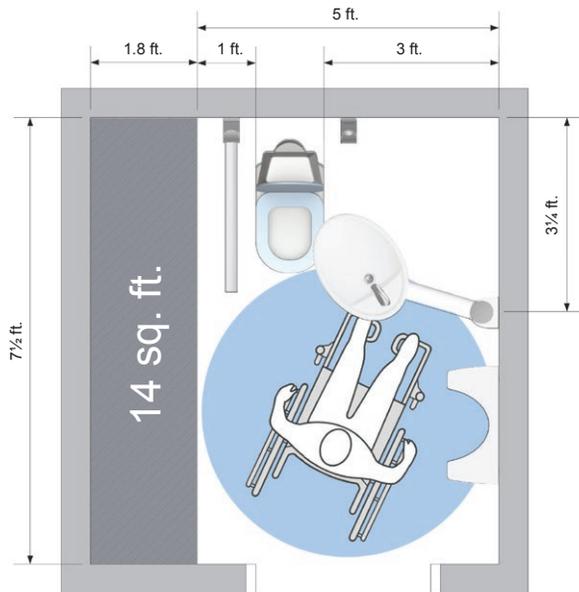


**Toilet for self-sufficient wheelchair users**  
With Ropox solutions the bathroom area can be reduced by 14 sq. ft. compared to a standard bathroom.

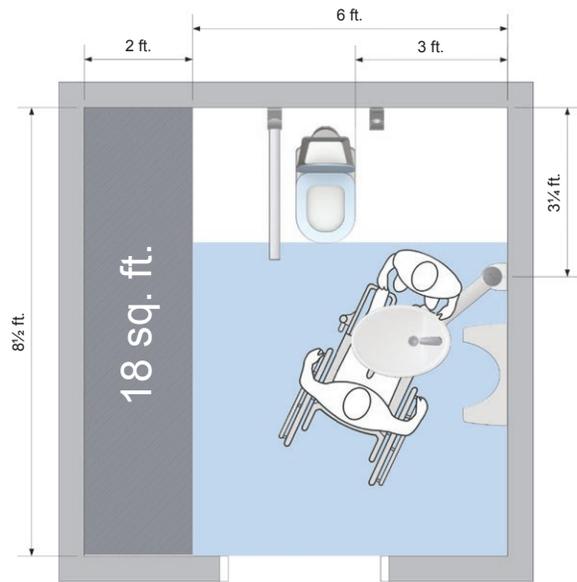

**Toilet for wheelchair users who need a caregiver's assistance**  
With Ropox solutions the bathroom area can be reduced by 18 sq. ft. compared to a standard toilet bathroom.


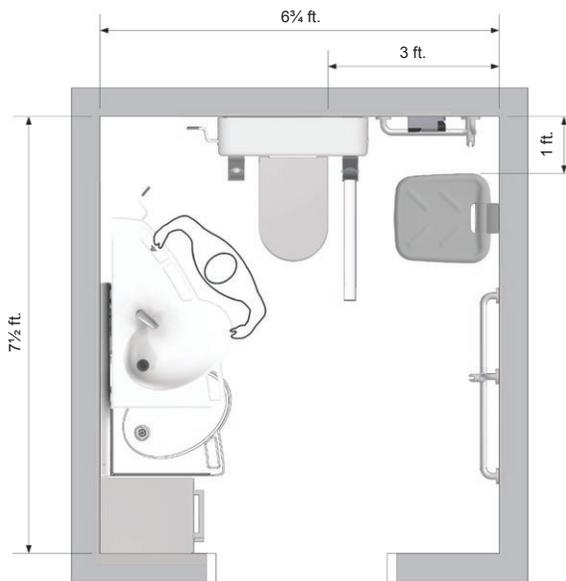


**Toilet for walking-impaired users**

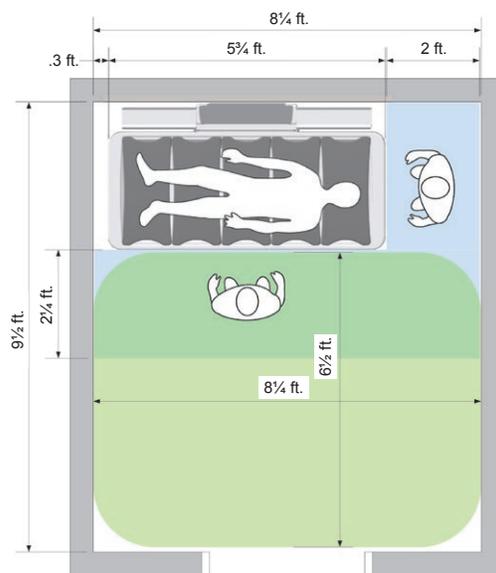

**Shower/changing bed for users who need a caregiver's assistance in institutions**

THE VARIANTS ARE AVAILABLE IN OUR DATA SHEETS  
WWW.PROCARE-INNOVATIONS.COM  
PHONE: 855-528-0421

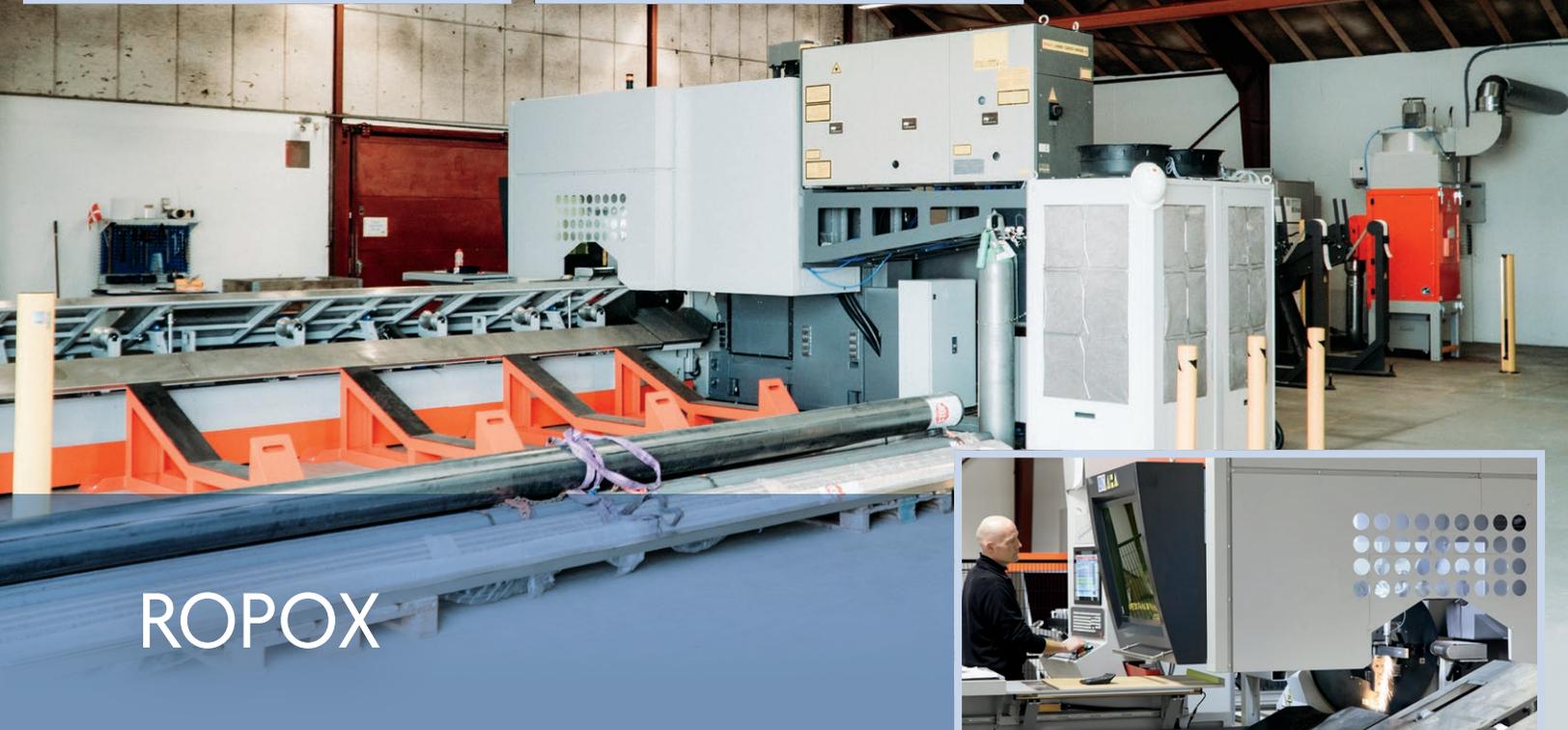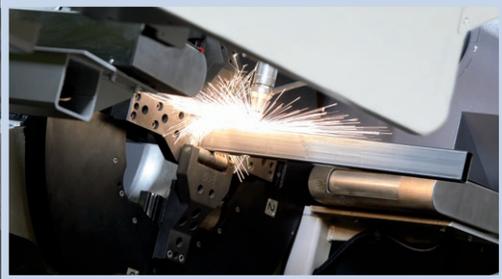

## BETTER WAYS TO BETTER DAYS

Ropox is a vertically integrated manufacturer of products, ranging from in-house design, to robotic welding and fabrication operations. With a focus on quality and detailed workmanship, the Ropox solutions have been creating Better Ways to Better Days for over 40 years.





BEAUTIFUL HEIGHT-  
ADJUSTABLE TABLES FOR  
INDIVIDUAL NEEDS



FLEXIBLE AND  
HEIGHT-ADJUSTABLE  
KITCHENS



INNOVATIVE AND  
EXCITING PRODUCTS  
IN A BEAUTIFUL  
COHESIVE SETTING



SCAN TO  
VIEW VIDEOS

759 Flory Mill Road, Lancaster PA 17601  
Toll Free: 855-528-0421 | [www.ProCare-innovations.com](http://www.ProCare-innovations.com)

