

”Ideal washbasin when space is limited”



ROPOX QUICKWASH

A HEIGHT-ADJUSTABLE WASHBASIN THAT CAN PROJECT OFF THE WALL

QuickWash is a height-adjustable washbasin with the ability to project off the wall. The sink can be lowered 14 cm (5.5”) from its fixed point and move 10 cm (4”) out from the wall.

It is an ideal solution for wheelchair users because it gives them the much-needed additional room for the legs. When the sink is projected out and down, wheelchair users get 63 cm (25”) legroom under the sink, allowing them to get close to the washbasin and reach the tap.

The height is adjusted with a release handle which also can be used as a towel holder. It requires very little effort to lower and raise the washbasin. The position lock secures that the basin is securely fastened.

QuickWash does not require electricity, and it comes preassembled, making it easy to install. The cost-effective solution is perfect for small bathrooms, and the clean, simple design makes it a good fit for private homes.



ROPOX™

BETTER WAYS TO BETTER DAYS


- A PART OF AddLife®

QuickWash Basin

Part No.	DESCRIPTION
40-42121	QuickWashbasin Complete UDI 57075810014R6
Inc. Washbasin, faucet & plumbing kit	



Range of motion

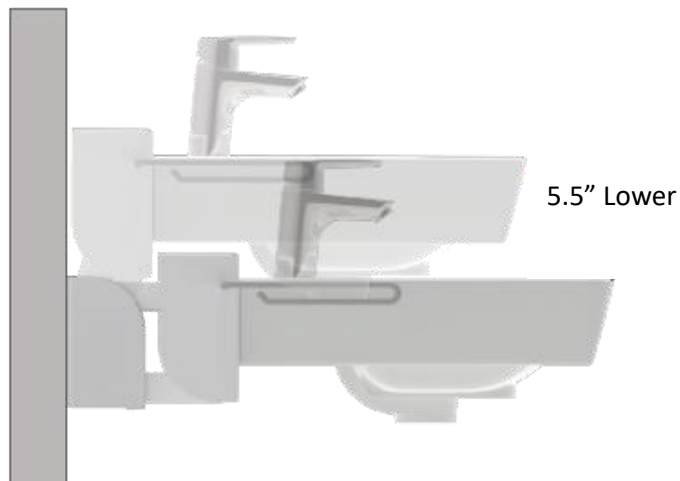

Technical data	
Max. user Weight	330 lb
Height adjustment	5.5 in
Forward movement	4.0 in
Complete weight	50.6 lb
MDR Class 2017/745	Class I
Expected Service Life	5 years (9000 Cycles) <i>Tested 18.000 cycles</i>
Intended environments	This product is to be used only in: Professional Healthcare Environment Home Care Environment