

5900 8200 Lotus / 5900 8250 Lotus compact

Lotus Active Standup Aid



TRANSFORMING CARE, IMPROVING LIFE

Lotus active standup aid with mechanical leg adjustment

The Lotus is an active standup aid designed for individuals who can still stand independently and possess sufficient core stability but may have difficulty walking or are unable to walk. Caregivers can safely and efficiently move the client after they have seated themselves on the Lotus.

The Lotus is intended for short transfers, such as moving from the bed to the toilet. Because the client assumes an active seated position, less effort is required to return to a standing position. Optimal comfort is provided for both the client and the caregiver. The Lotus features swiveling seat components, comfortable lower leg support, and base adjustment. This feature allows the Lotus to be brought close to the client, enabling them to actively take a seat. The Lotus encourages the client to actively participate in transfers, which helps train muscles and promotes independence.



5900 8200 Lotus

- Maximum load capacity of 440lb
- Rotating seat sections
- Mechanical Base Adjustment
- Wheel locking
- Variable hand grip positions
- Comfortable lower leg supports
- Ideal seating angle
- Sturdy footplate
- 2 lockable swivel wheels
- Adjustable leg support
- Antibacterial coating on frame.
- Smooth finish for optimal cleaning and disinfection.



Suitable for mobile people requiring some assistance.




5900 8250 Lotus compact



Scan to visit our website and learn more.




Rotating seat sections.

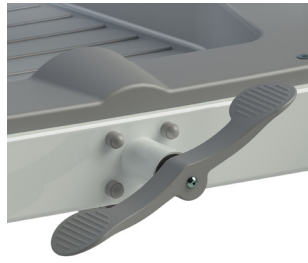


Mechanical base adjustment.



Variable contact points.adu

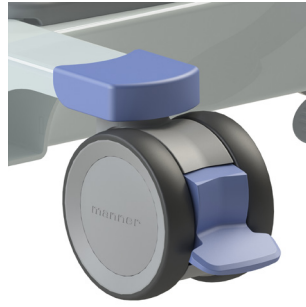


Wheel locking.

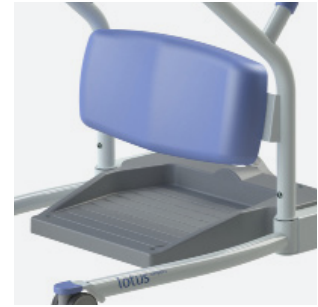


Comfortable lower leg support.

Technical specifications

Product no.	Product description	Material	Lifting capacity kg/lbs	Length x width x height	Width legs	Mechanical base spreading	Weight kg	Weight lbs	Client length
5900 8200	Lotus	steel	200 kg 440 lbs	96 cm x 68 cm x 110 cm 37,8" x 39,4" x 43,3"	68 cm 26,8"	68 cm - 96 cm 39,4" - 37,8"	32	77	160 cm - 195 cm 63" - 77"
5900 8250	Lotus compact	steel	200 kg 440 lbs	96 cm x 68 cm x 101 cm 37,8" x 39,4" x 39,8"	68 26,8"	68 cm - 96 cm 39,4" - 37,8"	32	77	130 cm - 165 cm 51" - 65"

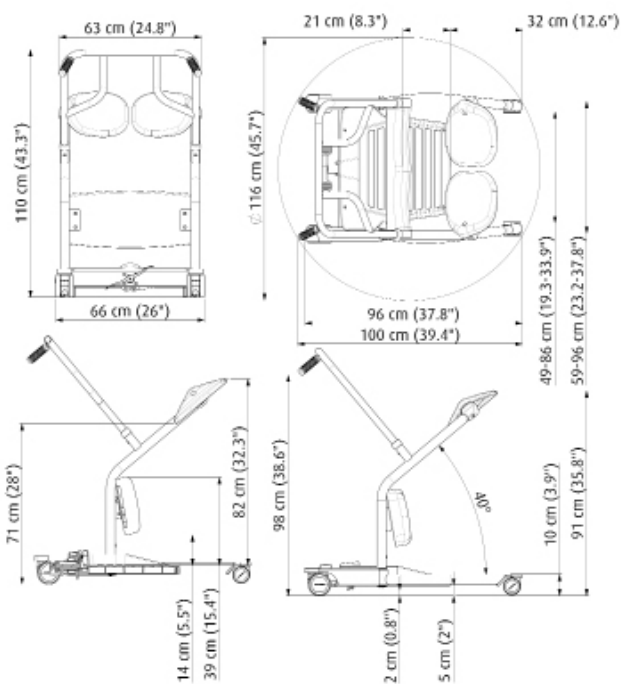

Lotus

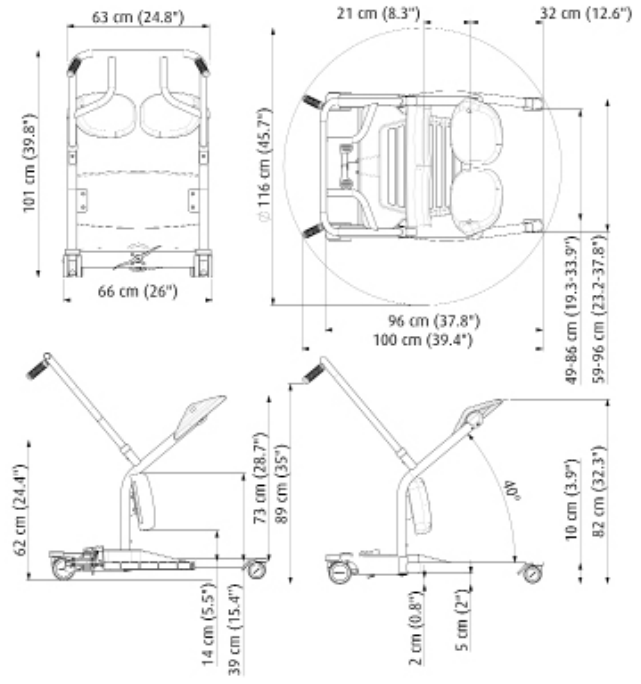

Lotus compact